

ClearVac

Disposable Smoke Evacuation Pencil

CLEAR VISION , SURGICAL PRECISION

★ **PATENTED TECHNOLOGY**

★ **ENHANCED VISIBILITY**

★ **SURGICAL SAFETY**

The **ClearVac Disposable Smoke Evacuation Pencil** delivers precise electrosurgical performance while efficiently capturing smoke at the source, ensuring a clear surgical view and optimal safety for the entire team.

Key Features

Patented Extendable Nozzle

Helical rotation, extends up to 20mm for precise smoke capture at any cutting depth.

Removable Transparent Suction Tip

Provides a clear surgical view and enables quick electrode replacement.

Patented Buffer Ball Joint

The inertia buffer absorbs 85% of sudden forces.

360° Swivel Connector


Allows effortless maneuverability.

Water-Resistant PCB

Sealed PCB design prevents liquid ingress for enhanced durability and safety.

Ergonomic Grip

Lightweight ribbed handle minimizes fatigue and enhances control.

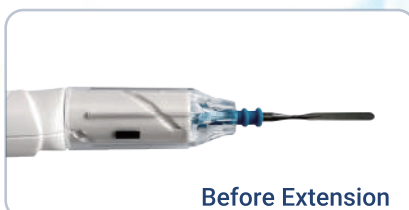

Promed Supply Company Ltd

Flat A, 10/F, Mai Wah Ind. Bldg., 1-7 Wah Sing St, Kwai Chung, N.T, Hong Kong



www.the-promed.com




info@the-promed.com




+852 9385 0330

Clinical Benefits

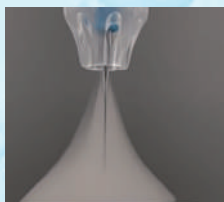


Optimal Smoke Capture

- **Extendable Nozzle** extends up to 20 mm for consistent smoke capture at any surgical depth.
- **Enlarged Evacuation Chamber** with optimized suction holes enables fast, efficient smoke removal.
- **Transparent Tip** ensures clear visibility of the surgical site.

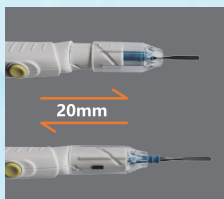


Maximum Safety Engineered

- **Patented Inertia-Buffer Joint** absorbs sudden movements for precise control.
- **Sealed Waterproof PCB** prevents liquid ingress and eliminates electrical failure risk.
- **X-Ray Detectable Tip** enhances patient safety.

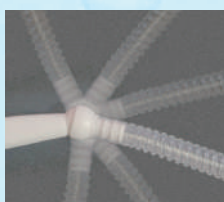

Ergonomic Precise Control

- **Ergonomic Ribbed Grip** delivers superior control and tactile feedback for a secure hold.
- **360° Swivel Connector** enables smooth, unrestricted movement for precise surgical performance.
- **Universal Compatibility** with standard electrosurgical units and hex electrodes.

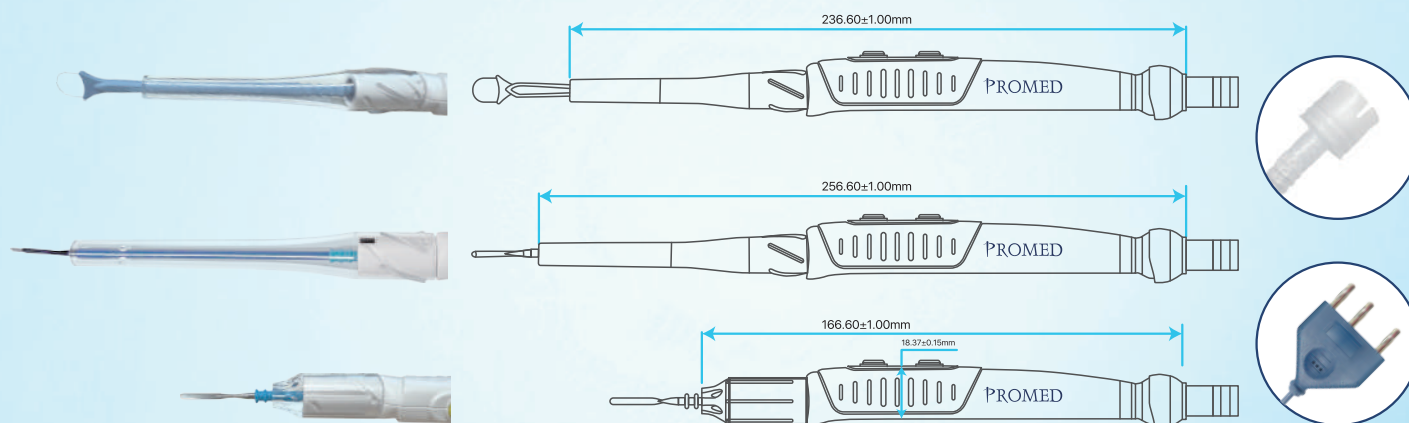

Specifications

Product Code	Electrode	Hose Ø	Hose Length	Cable Length	Plug	Packing
PD900B	Ø2.36 69mm Stainless Steel Blade	22mm	Eva Tube 3M	3M	3 Pin Plug	20pcs/Case
PD900B5	Ø2.36 69mm Stainless Steel Blade	22mm	Eva Tube 4.5M	4.5M	3 Pin Plug	20pcs/Case
PD900BT	Ø2.36 69mm Non-Stick Blade	22mm	Eva Tube 3M	3M	3 Pin Plug	20pcs/Case
PD900B5T	Ø2.36 69mm Non-Stick Blade	22mm	Eva Tube 4.5M	4.5M	3 Pin Plug	20pcs/Case
PD120R	Ø2.36 EL2V65 Loop Blade	22mm	Eva Tube 3M	3M	3 Pin Plug	20pcs/Case
PD120R5	Ø2.36 EL2V65 Loop Blade	22mm	Eva Tube 4.5M	4.5M	3 Pin Plug	20pcs/Case
PD110B	Ø2.36 152mm Stainless Steel Blade	22mm	Eva Tube 3M	3M	3 Pin Plug	20pcs/Case
PD110B5	Ø2.36 152mm Stainless Steel Blade	22mm	Eva Tube 4.5M	4.5M	3 Pin Plug	20pcs/Case
PD110BT	Ø2.36 152mm Non-Stick Blade	22mm	Eva Tube 3M	3M	3 Pin Plug	20pcs/Case
PD110B5T	Ø2.36 152mm Non-Stick Blad	22mm	Eva Tube 4.5M	4.5M	3 Pin Plug	20pcs/Case



Promed Supply Company Ltd

Flat A, 10/F, Mai Wah Ind. Bldg., 1-7 Wah Sing St, Kwai Chung, N.T, Hong Kong



www.the-promed.com



info@the-promed.com



+852 9385 0330